

VBM Endoscopy Mask

Allows simultaneous ventilation and endoscopy of the patient to increase patient comfort and the working quality of the doctor.

Application:

- **Fiberoptic Intubation**
The ET Tube can be railroaded orally or nasally along the bronchoscope through the flexible membrane. After intubation the membrane is pushed inside which allows easy removal of the mask.
- **Bronchoscopy**
The membrane guarantees an ideal seal and therefore an efficient ventilation before, during and after bronchoscopy without time pressure.
- **Gastroenterology**
- **Trans Esophageal Echocardiography**

The use of Entonox (50% N₂O, 50% O₂) is ideal for minor procedures where pain is predictable and of short duration. Entonox can be self administered from the conscious patient via the mask. It has a rapid analgesic effect and therefore allows a pain free insertion of the Endoscope through the hole in the membrane. Entonox wears off quickly without post operative problems.



Fiberoptic Intubation




Endoscopy Mask

Fiberoptic Intubation / Bronchoscopy

Mask	Patient	Hole	REF	Box
# 1	Infant	2mm	30-40-111	6
# 3	Child	3mm	30-40-333	6
# 5	Adult	5mm	30-40-555	6

Gastroenterology / Trans Esophageal Echocardiography

Mask	Patient	Hole	REF	Box
# 5	Adult	10mm	30-40-777	6

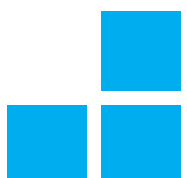

Fixation Tapes



Harness

Material: silicone

Size	REF	Box
Child	35-70-155	10
Adult	35-70-255	10

Hook Ring

for attachment between mask and membrane

REF 80-11-266 Box 10



Fixation Tape

Material: foam

For Mask #	REF	Box
# 1	35-40-622	10
# 3+5	35-40-633	10

Membrane

Silicone Membrane

autoclavable up to 134°C

Hole	REF	Box
2mm	30-40-101	5
3mm	30-40-103	5
5mm	30-40-105	5
10mm	30-40-110	5



Bronchoscope Airway

The VBM Bronchoscope Airway avoids damage to the bronchoscope by the patient.

1. The lumen is large enough to guide the bronchoscope and Endotracheal Tube
2. The length allows to guide the bronchoscope into the pharynx
3. The external biteblock protects the bronchoscope
4. For single use – no risk of infection – no risk of contamination
5. Different sizes allow adjustment to patient and application

Bronchoscope Airway

with external biteblock

single use, Material: EVA

#	Patient	Length	REF	Box
0	Newborn	36mm	30-40-400-1	10
2	Child	65mm	30-40-420-1	10
4	Adult	85mm	30-40-440-1	10

